


tg barbed broaches

SKU: 1640

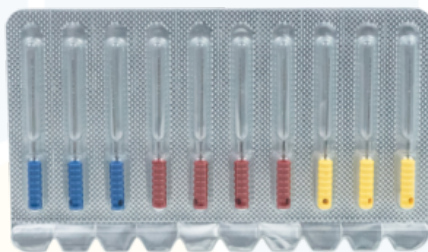

CE 0123

Root Canal Preparation

tgbarbed Broaches

- Exceeds all ISO/ADA specifications
- High Quality Stainless Steel
- Hand use
- Blister card packing 10 pieces/blister in working lengths of 21mm
- Color coded, ISO colors & sizes
- CC-Cord Handle
- Sterile
- About 40 Flexible barbs are spirally arranged on the coned shaped core

- The size of the barbs corresponds to half of the core diameter
- The delicate barbs are resilient and will not break in the canal if used properly
- Sterilizable up to 200 C /390 F

Available Sizes

- XXXXF (410), XXXF (420), XXF (430), XF (440), F (450), M (460), C (470), assort. XXF-F (480)
- Length 21mm



Sales & Information

Technical & General, 2 Albion Place, London W6 0QT United Kingdom
Tel: +44 20 8834 7575 Email: sales@technical-general.com
www.tgdent.com

Technical & General Ltd - Generated on 11.09.26